

## Fiber-optics for Optogenetics

For *in-vivo* and *in-vitro* stimulation

Ver. 03

### Introduction

Prizmatix provides a full range solution to the optogenetics *in vivo* and *in vitro* fiber optics. Made of silica or polymer fibers, the high NA fibers are assembled to fit any research set-up with various combinations of connectors, ferrules, core diameters and lengths.



### *In-vivo* fibers - example


This is a Y-shaped fiber for simultaneous *in-vivo* stimulation of two hemispheres in behaving animals.

Two high NA flexible polymer fibers are accommodated in an FC/PC connector, connected to the output port of a rotary joint. The two fiber branches are terminated with ferrules, to connect to the fiber cannulae implanted in the animal's head.



### *In-vitro* fibers - example

This is a single fiber with a thin stainless steel tip at its distal end. A high NA 200 $\mu$ m silica fiber is protruding from the tip for several millimeters. The stainless steel enables anchoring the fiber to a micromanipulator for precise movement and spatial location control. Yet, since the outer diameter of the stainless steel is 3-4 times the core diameter of the fiber, it would have been obscured the image under the microscope. Therefore, we design and build these fibers with the polished bare fiber protruding from the tip for several millimeters, improving the positioning accuracy of the fiber over the target in the brain slice.

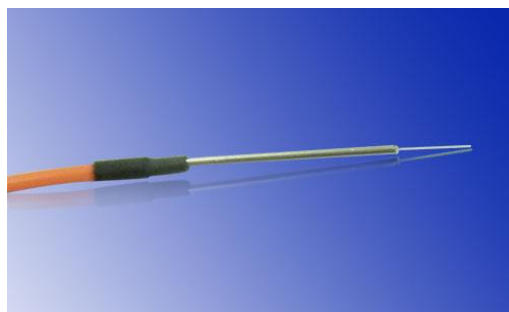

#### Main Office

Phone: +972-72-2500097  
Fax: +972-72-2500096  
[sales@prizmatix.com](mailto:sales@prizmatix.com)

#### European Sales Office

Phone: +44-(0)77-9172-9592  
Fax: +44-(0)20-7681-2977  
[sales.europe@prizmatix.com](mailto:sales.europe@prizmatix.com)

#### North America Sales Office

Phone: +1-(248)-436-8085  
Fax: +1-(248)-281-5236  
[sales.usa@prizmatix.com](mailto:sales.usa@prizmatix.com)

P.O.B. 244 Givat-Shmuel 5410102, Israel

# Prizmatix

## Optional Parameters

- **Fiber specification:**
  - Polymer optical fibers. NA 0.5. Core diameters [ $\mu\text{m}$ ]: 250, 500, 750, 1000, 1500
  - Polymer optical fibers. NA 0.63. Core diameters [ $\mu\text{m}$ ]: 500, 1000, 1500
  - Silica fibers NA 0.53. Core diameters [ $\mu\text{m}$ ]: 200, 400
  - Silica fibers NA 0.37. Core diameters [ $\mu\text{m}$ ]: 100
- **Fiber length:** as required, usually 1-1.5 meters.
- **Connectors:** SMA, FC/PC
- **Number of branches:** 1, 2, 3.
- **Ferrules material:** Stainless, Zirconia
- **Ferrules core [mm]:** 2.5, 1.25, Custom
- **Stainless steel tip diameter [mm]:** 0.5 - 2.0
- **Stainless steel tip length [mm]:** 10 - 60
- **Fiber protruding length [mm]:** 1 - 10

**SPECIFY PARAMETERS WHILE ORDERING**

## Optional Accessories

### Ultra High Power LEDs (UHP LEDs):

The Ultra High Power LED is a light source for fluorescence microscopy and other applications, and is an effective replacement of spectral lamps and lasers in many applications. The LED driver supports CW or external TTL modulation with user-controlled frequency and duty cycle. For more details please see our website:

<http://www.prizmatix.com/UHP/ultra-high-power-leds.htm>

### Beam Switcher:

The Beam Switcher accessory allows for Prizmatix LEDs installed on a microscope to be used either as microscope epi-fluorescence illumination or illuminating a specimen via a fiberoptic probe.

Especially useful for in-situ Optogenetics studies, the Beam Switcher saves cost, space and allows maximum flexibility in experiment design. For more details refer to

Prizmatix Website: <http://www.prizmatix.com/optogenetics/beam-switcher.htm>

### Fiber Coupler and Liquid Light Guide Adaptors:

The side port of the Beam Switcher can be easily configured to fiber coupling by means of the FCA Fiber Coupler Adaptor (SMA, CF or ST connector). This adaptor can be easily assembled by a user on the Mic-LED output.

Please see video clip <http://youtu.be/iv7dlwLHaUE>.

---

#### Main Office

Phone: +972-72-2500097  
Fax: +972-72-2500096  
sales@prizmatix.com

#### European Sales Office

Phone: +44-(0)77-9172-9592  
Fax: +44-(0)20-7681-2977  
sales.europe@prizmatix.com

#### North America Sales Office

Phone: +1-(248)-436-8085  
Fax: +1-(248)-281-5236  
sales.usa@prizmatix.com

---

P.O.B. 244 Givat-Shmuel 5410102, Israel

# Prizmatix

---

## Filter Wheel:

LEDs can be equipped with a 6-positions filter wheel at the beam output, especially useful for the UHP-Mic-LED-White light source. For more details please see our website: <http://www.prizmatix.com/uhp/uhp-mic-led-white.htm>


## Rotary Joint:

Prizmatix Rotary Joint is especially designed for Optogenetics experiments with High NA fibers equipped with FC connectors. Rotary Joint very low friction and smooth swivel operation are demonstrated in this video -

<http://www.youtube.com/watch?v=q5RE14fove0>, showing continuous rotation for over 30 seconds after a fingertip startup.

The torque required for free movement is very low thus reducing the Rotary Joint's behavioral effect on free moving animals. Prizmatix Rotary Joint can be used with large core input optical fiber connected to several thin output fibers simultaneously. This allows concurrent delivery of light to separate areas of the brain without loss of power due to the fiber split. For more details refer to Prizmatix Website:

<http://www.prizmatix.com/optogenetics/rotary-joint.htm>

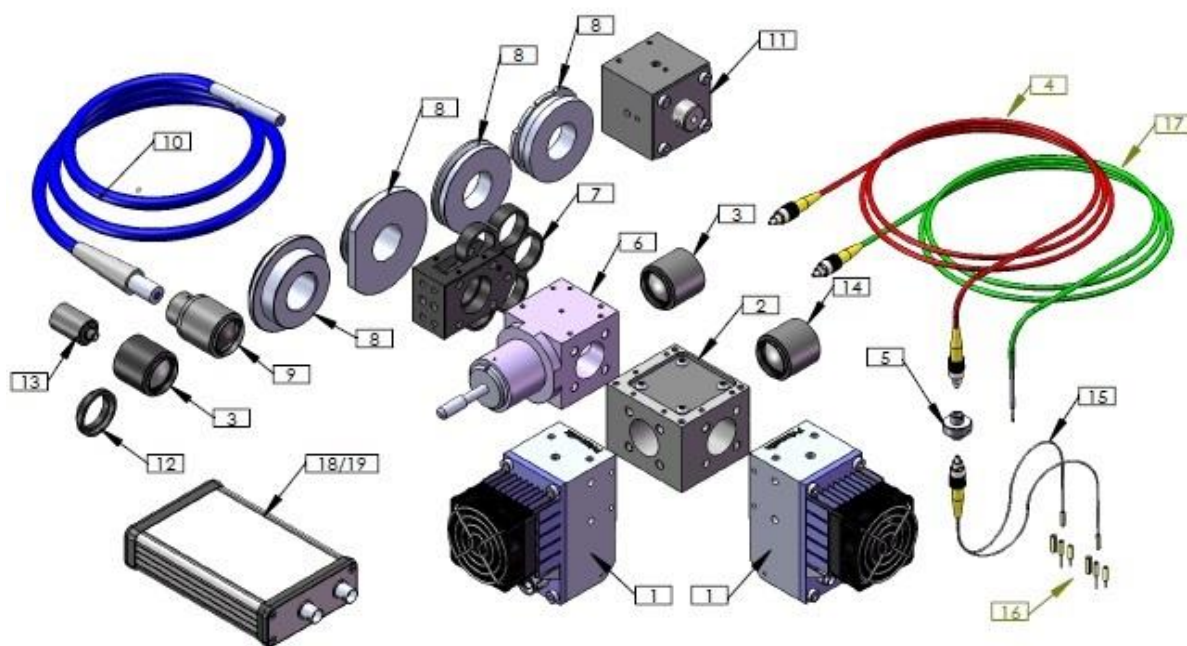

---

### Main Office

Phone: +972-72-2500097

Fax: +972-72-2500096

[sales@prizmatix.com](mailto:sales@prizmatix.com)

### European Sales Office

Phone: +44-(0)77-9172-9592

Fax: +44-(0)20-7681-2977

[sales.europe@prizmatix.com](mailto:sales.europe@prizmatix.com)

### North America Sales Office

Phone: +1-(248)-436-8085

Fax: +1-(248)-281-5236

[sales.usa@prizmatix.com](mailto:sales.usa@prizmatix.com)

---

P.O.B. 244 Givat-Shmuel 5410102, Israel